





The Clinic-Chek has been designed for simple analysis of specific body fluids in medicine. Three scales provide measurement of serum protein levels and urine specific gravity, also a refractive index scale is provided when a cross reference to a specific fluid chart is called for.


SIMPLE SET UP:

1. Insert the included batteries with positive side facing outward.

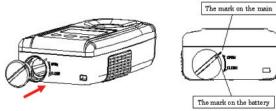

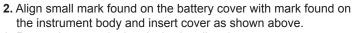

3. Rotate the cover in the clock wise direction until it is in the position shown to the right.

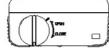

CLEANING:

Cleaning of the measurement surface and well should be performed immediately after each sample reading. Never immerse the instrument in any liquid. When the measurement surface and well have been completely cleaned no residue should be present.

To properly clean the sample well and glass measurement surface use a mild soap and water solution or Isopropyl Alcohol followed by a distilled water rinse and then thoroughly dry with a soft lint and residue free cloth or a product such as Kimwipes®.

To clean the refractometer's body use a soft cloth with a mild cleaning solution like window cleaner applied to the cloth and follow by wiping the instrument dry.

CLEANING PRECAUTIONS:

- 1. Never use any harsh cleaning agents that will damage the instrument.
- 2. Strong solvents should never be used. Such chemicals will deteriorate the prism seal and attack individual component material which will adversely affect or destroy the instrument.

CALIBRATION:

- 1. Be sure that the measuring surface and well are clean as outlined under the "Cleaning" section.
- 2. Apply distilled water to the measurement surface.
- 3. Allow time for the distilled water sample to temperature equilibrate to that of the instrument.
- 4. Press and hold the "CAL" key until CAL is displayed. While still holding the "CAL" key press the "READ" key and then release both kevs.
- 5. A series of dashed lines will appear and successful calibration will be indicated when the word END is displayed.

SAMPLE ANALYSIS:

Before analyzing a sample, it should be noted that to attain the greatest accuracy the refractometer should be calibrated at the same temperature as the environment that it will be used in. As an example; if the environment that the instrument is being used in changes by five degrees, then the instrument should be re-calibrated.

- 1. Inspect the glass measuring surface and well to be sure that there is no residue remaining from a previous sample analysis. If the measuring surface and well does not appear to be clean, then refer to the cleaning section of this user guide before proceeding.
- 2. Extract the sample to be tested and apply adequate sample to the sample well to completely cover the glass measuring surface. Minimum sample size is 0.3 mL.
- 3. Allow the sample time to reach the same temperature as the instrument. Failure to do so will result in inaccurate measurement.
- 4. Press the "READ" key and note the displayed value and the position of the annuciator. The annunciator should be located below the fluid scale that you are attempting to attain a reading for. If the position of the annunciator is not located under the desired scale, then press the "SCALE" key until the annuciator position indicates the correct scale.

SCALES:

To change the reading mode scale between Urine specific gravity. Serum protein or refractive index (nD) press the "SCALE" key to toggle through those choices. The annuciator displayed on the screen below the scale function denotes the active scale.

Displaying Actual Temperature:

To display the temperature at the measuring surface, momentarily press the "CAL" button while the sample value is displayed.

AUTOMATIC TEMPERATURE COMPENSATION:


All Reichert refractometers are automatically temperature compensated assuring extremely accurate results. Automatic temperature compensation corrects readings over a range of temperatures. As an example; samples taken within the working temperature range of the instrument are corrected to a reference temperature of 68° F (20° C). Temperature compensation is essential because refractive index varies inversely with temperature. This feature is used in lieu of temperature control of the measuring surface and sample. For this feature to be effective however, it is necessary to allow the sample to temperature equilibrate to the ambient temperature of the measuring surface. Without allowing for temperature equilibration of the sample the read values will not be as accurate as they could be since the sample will be in a state of temperature and refractive index change.

ERROR CODES:

- Err01 No sample present. Add sample.
- Err02 Inadequate sample. Add additional sample.
- Err03 Sample exceeds the refractive index reading range.
- Err04 Sample temperature has not equilibrated. Allow more time for temperature equilibration.
- Err05 Excessive ambient light. Cover sample well when reading.
- Err06 Excessive ambient light. Cover sample well when reading.
- Err07 Positive calibration error. Re-calibrate with distilled water.
- Err08 Negative calibration error. Re-calibrate with distilled water.
- Err09 Poor sample condition. Sample may not be able to be read.
- Err10 Sample type may not be able to be read.
- Err12 Index of sample is out of range of instrument.
- Err5X For any errors in this series contact Reichert Analytical Instruments for technical assistance.

BATTERY CONDITION:

The following display examples represent either a low battery condition or an extremely low battery condition where the instrument will no longer operate properly. In either case the batteries should be promptly replaced.

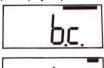

Battery condition low. Replace Batteries.



Battery condition extremely low. Replace Batteries.

SPECIFICATIONS:

Catalog 13940021

Reading Scales - Urine specific gravity, USG Serum Protein, g/100mL Refractive index, nD

Reading Ranges,

USG, 1.000 to 1.040, Accuracy +/- 0.005 Serum Protein, 1.6 to 15.0, Accuracy +/- 0.1 Refractive index, 1.3330 - 1.4465, Accuracy +/- 0.0002

Calibration - Distilled Water

Automatic Temperature Compensation - 68°F (20°C)

The built in temperature compensation mode is based on the temperature coefficients that have been calculated for the individual solutions offered in this instrument.

Prism - Glass

Illumination - 589nm LED

Dimensions - 54 x 27 x 100 mm / 2.13 x 1.06 x 3.9 inches

Power - 2 AAA Batteries, included.

Battery life - 10,000 readings, Auto Off Sleep Mode.

- Ratings IP65 Dustproof/Water Resistant, CE, RoHS, and WEEE compliant.
- Warranty One year against manufacturing defects. Evidence of tampering voids warranty.



Reichert Technologies Analytical Instruments 3362 Walden Ave Depew, NY 14043 USA

Toll Free: 888-849-8955 Phone: 716-686-4500 email: reichertai.refractometers@ametek.com www.reichertai.com

> AMETEK GmbH Business Unit Reichert Carl-von-Linde-Strasse 42 85716 Unterschleissheim Munich, Germany Tel: +49 (89) 315 8911 0 Fax: +49 (89) 315 891 99

Œ

ISO-9001 Certified